Los Alamos National Laboratory is operated by the University of California for the United States Department of Energy under contract W-7405-ENG-36

LA-UR--86-3119

DE87 000142

RIMS DIAGNOSTICS FOR LASER DESORPTION/LASER ABLATION TITLE:

E. C. Apel, N. S. Nogar, C. M. Miller, and R. C. Estler AUTHOR(S):

Presentation at the Third International Symposium on SUBMITTED TO: Resonance Ionization Spec. oscopy & Its Applications, September 7-12, 1986, Swansea, Wales, UNITED KINGDOM

(To be published in Conference Proceedings)

DISCLAIMER

This report was prepared as an account of work sponsored by an agency of the United States Consermment Neither the Pasted States Government nor any agency thereaf, nor any of their employees, makes any warranty, express or mighted, or assumes any legal limbility or responsibility to the accuracy, completeness, or usefulness of my information, apparatus, product, or process disclosed, or represents that its use would not infringe privately owned rights. Refer ence betein to any specific commercial product, process, or service by trade name, trademark, manufacturer, or otherwise does not necessarily constitute or imply its endorsement, recommendativa, or favoring by the United States Government or any agency thereof. The views and anomous of authors expressed begain the not necessarily state or reflect those of the United States Government or any agency thereof

By acceptance of this accept th - a nonexclusive, royally-live license to publish or reprinduce the published from of this combitment on to allow others to do so, for U.S. Government purposas

The Los Almins National Laboratory requests that the publisher identify this article as work performed under the auspices of the U.S. Department of Livergy

S Alamos National Laboratory Los Alamos, New Mexico 87545 (111)



1 131M NO. 8 M 114 51 No. 2829 5-81

DELY REDUCTION OF THIS DOCUMENT IS UNLIMITED

RIMS Diagnostics for Laser Desorption, Laser Ablation

E. C. Apel, N. S. Nogar, and C. M. Miller and R. C. Estler'

Los Alamos National Laboratory, Los Alamos, New Mexico 87545, 'Ft. Lewis College, Durango, Colorado 81301

Abstract. Laser desorption mass spectrometry is a useful method for interrogating materials and events at or near surfaces. Laser desorption ablation combined with Resonance Ionization Mass Spectrometry (RIMS) provides a powerful tool to obtain information on chemical composition and speciation and, in some cases, internal and translational energy distributions. The application of this technique to the interrogation of materials and interfaces is discussed for several systems, including the analysis of conventional analytical samples, and the study of optical damage events.

1. Introduction

RIMS is a powerful tool for elemental and isotopic analysis. The multistep laser photoionization process, when coupled with conventional mass analysis, can provide exceptional performance in detectivity, dynamic range, and discrimination against interfering species (Nogar et al 1985a, Fassett et al 1985). These properties can be used to great advantage both in the analysis of conventional materials (Bekov and Letokhov 1983) and in the interrogation of interfacial phenomena (Kimock et al 1983). In this report we will discuss the use of laser desorption ablation in combination with RIMS analysis.

The motivations for this work are threefold. First, we desire to improve the effective duty cycle, relative to continuous sample evaporation, for the examination of routine analytical samples. Second, this method will allow the direct analysis of materials without extensive sample preparation. And third, we wish to study the mechanism of laser-material interactions, particularly optical damage, and characterize the desorption, ablation/damage process.

2. Experimental

The apparatus has been described previously (Nogar et al 1985b, Estler et al 1986), and will only be summarized here. Briefly, a Q-switched Nd:YAG laser equipped with beam-filling optics is used for desorption of a sample mounted on a manipulator in the source region of a time-of-flight mass spectrometer. Typical

laser parameters include 10 nsec pulse length, 5 to 50 mJ pulses, 0.2 to 1 mm spot size at the sample, and 10 llz repetition rate. After a variable time delay, the interrogation laser (XeCl excimer-pumped dye laser, typical parameters: ~12 nsec pulse, I mm diameter, 0.5 to 3 mJ. ~0.3 cm⁻¹ bandwidth) is passed through the spalled plume. At a fixed time delay following the excimer trigger, ions are detected with a channel electron multiplier, and the signal processed with a gated integrator and standard analog signal processing electronics.

3. Results and Discussion

Two-photon excitation. In order to expedite diagnostics for our laser desorption ablation studies, we wished to minimize the complexity of the laser ionization process and associated hardware. We have therefore used single-color ionization processes at the fundamental output frequency of the interrogation dye laser. This, in turn, required n-photon excitation, n=2, for moderate to high ionization potential elements. We have examined (Apel et al 1986) a number of atomic (and some molecular) systems ionizable via "2-1" processes. A list of observed transitions is shown in Table 1. For the atomic systems, the two-photon excitation typically saturated at ~2 mJ pulse, while the ionization step required somewhat higher energies.

TABLE 1
Observed Two-Photon Transitions

	Ground State		Excited State		Energy Difference
Element	Configuration	Term	Configuration	Term	(cm ⁻¹)
Ca	$(4s^2)$	18	(4s4d)	1 D ₂	37298
CaF		$X_{-2}\Sigma$.		$\mathbf{F}^{-2}\pi$	37550
Ta	(5d ³ 6s²)	4F3 2	(5d ³ 6s7s)	J 3.'2	44096
	$(5d^36s^2)$	4Fa 2	$(5d^{3}6s7s)$	${}^4{ m F}_{3,2}$	43964

^aT., is given for CaF.

Analytical Applications. The use of pulsed lasers with continuous sample evaporation can result in substantial loss of analyte because of the low effective duty cycle. For a probe laser beam diameter of 0.5 cm, and an atomic velocity of 5 x 10⁴ cm sec, the rate of sample turnover in the beam volume is $\sim 10^5$ sec. ¹. For a laser repetition rate of 10 sec. ¹, the effective duty cycle is $\sim 10^{-4}$. The use of lasers (or particle beams, Kimock et al 1984) to pulse desorb the sample can substantially improve this value. We have recently (Nogar et al 1985b) demonstrated laser desorption RIMS for tantalum samples. Tantalum atoms were detected via a "2 · 1" ionization through the "F₃₋₂ state. The results were quite promising: the pulse of desorbed material was sufficiently narrow (· 4 μ sec) that the effective duty cycle was improved to $\sim 10^{-1}$, resulting in a substantial increase in sample utilization efficiency. The observed energy distributions were somewhat

anomalous: for high intensities (10° W cm²) the hydrodynamic temperature was measured to be 8000 K, while the kinetic temperature was 400 K, and the internal (electronic) temperature was ~2000 K. At lower intensities (4 x 10° W cm²) the observed distributions were cooler, and more nearly thermal, as shown in Fig. I. These results are consistent with the existance of a thermal desorption barrier (Nogar et al 1985b).

Laser Damage Studies. RIMS has also been used as a monitor of the interaction between lasers and dielectric materials (Estler et al 1986). In initial experiments, optical damage on uncoated CaF₂ substrates was initiated with 1.06 µm pulses at fluences of I to 10 J cm². Interrogation of the spalled plume revealed Ca atoms and CaF radicals only when damage events occured. Both species are monitored by "2 - 1" ionization processes through the ${}^{1}\mathrm{D}_{2}$ and F ${}^{2}\pi$ states. respectively. This choice of intermediates states was made so that both species could be detected with a single laser dye. In subsequent experiments, CaF₂ damage v as also induced at 355 nm and 266 nm. For 1.06 μ m irradiation, we observed thermal (850 K) velocity distributions for both Ca and CaF. In addition. the CaF radical exhibited significant amounts of internal

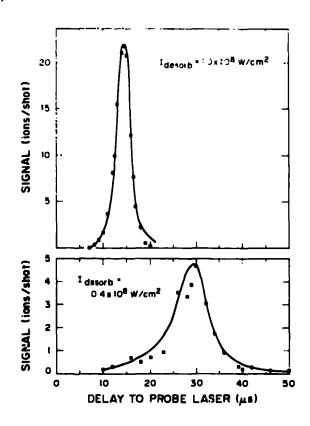


Fig. 1. Tantalum atom flight time distributions for high intensity (~10⁸ W cm²) and low intensity (~4 x 10⁷ W cm² laser desorption.

(rotational and translational) excitation. For both 355 nm and 266 nm irradiation, the velocity distributions were bimodal, with a fraction ($^{\circ}50\%$) of the spalled material exhibiting very high (4000 K) kinetic temperatures (see Fig. 2), while the remainder exhibited a temperature similar (800 to 1000 K) to that observed for the 1.06 μ m experiments. In addition, both the vibrational and rotational temperatures of the CaF radicals decreased with decreasing damage wavelength. Lastly, the threshold for damage decreased slightly with decreasing wavelength.

These results suggest that while the 1.06 μ m experiments can be adequately modeled in terms of a single damage mechanism (likely avalanche breakdown), the short-wavelength results suggest the onset of second mechanism, perhaps multiphoton absorption. This is consistent with both the bimodal velocity distribution, and with the decrease in CaF vibrational and rotational excitation.

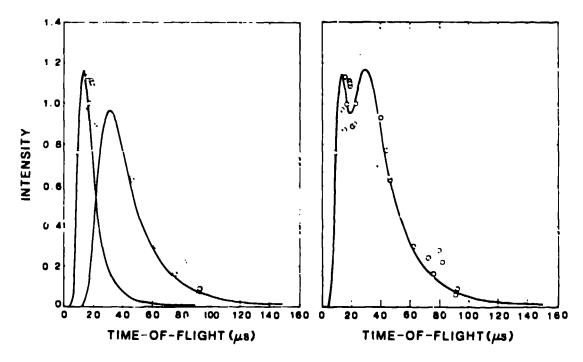


Fig. 2. (a) Calcium atom flight time distributions (dots) for damage induced at 266 nm. 25 J cm². The solid lines are Boltzmann fits for 4000 K and 850 K distributions. (b) Same as (a), where the solid line now represents a composite of the fast (33%) and slow (67%) contributions.

References

Apel E C. Anderson J E. Estler R C. Nogar N S. and Miller C M 1986 Appl. Opt. (submitted)

Bekov G I, and Letokhov V S 1983 Appl. Phys. B 30 161

Estler R.C. Apel F.C. and Nogar N.S. 1986 J. Opt. Soc. Amer. (accepted)

Fassett J D, Moore L J, Travis J C, and DeVoe J R 1985 Science 230 262

Kimock F M, Baxter J P, Pappas D L, Kobrin P II, and Winograd N 1984 Anal. Chem. 56 2782

Kimock F M, Baxter J P, and Winograd N 1983 Surf. Sci. 124 L41

Nogar N S. Downey S W, and Miller C M 1985a Spectroscopy 1 56

Nogar N S. Estler R C, and Miller C M 1985b Anal. Chem. 57 2441