

LA-9172

C.3

Los Alamos National Laboratory is operated by the University of California for the United States Department of Energy under contract W-7405-ENG-36.

CIC-14 REPORT COLLECTION
REPRODUCTION
COPY

LOS ALAMOS NATL. LAB. LIBS.
3 9338 00312 3030

*Spectrophotometric Determination of
Plutonium with
Chlorophosphonazo III in n-Pentanol*

Los Alamos Los Alamos National Laboratory
Los Alamos, New Mexico 87545

This work was supported by the US Department of Energy, Office of Safeguards and Security.

Edited by Arthur Nichols
Photocomposition by Pamela H. Mayne

DISCLAIMER

This report was prepared as an account of work sponsored by an agency of the United States Government. Neither the United States Government nor any agency thereof, nor any of their employees, makes any warranty, express or implied, or assumes any legal liability or responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product, or process disclosed, or represents that its use would not infringe privately owned rights. References herein to any specific commercial product, process, or service by trade name, trademark, manufacturer, or otherwise, does not necessarily constitute or imply its endorsement, recommendation, or favoring by the United States Government or any agency thereof. The views and opinions of authors expressed herein do not necessarily state or reflect those of the United States Government or any agency thereof.


LA-9172

UC-4

Issued: March 1982

Spectrophotometric Determination of Plutonium with Chlorophosphonazo III in n-Pentanol

N. M. Saponara
S. F. Marsh

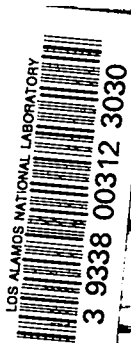

Los Alamos Los Alamos National Laboratory
Los Alamos, New Mexico 87545

SPECTROPHOTOMETRIC DETERMINATION OF PLUTONIUM WITH CHLOROPHOSPHONAZO III IN *n*-PENTANOL

by

N. M. Saponara and S. F. Marsh

ABSTRACT

Microgram amounts of plutonium are measured spectrophotometrically as the plutonium-chlorophosphonazo III complex after extraction into *n*-pentanol from 1.5*M* HCl. The relative standard deviation is 1.5% for the range of 2.5 to 17.5 μ g. The tolerance is excellent for many metals and nonmetals present in nuclear fuel-cycle materials. A preceding anion-exchange-column separation increases tolerance for certain metals and nonmetals.

I. INTRODUCTION

Chlorophosphonazo III (3,6-bis-[(4-chloro-2-phosphonophenyl)azo]-4,5-dihydroxy-2,7-naphthalenedisulfonic acid), termed CPA, has been used for spectrophotometric determinations of many elements including thorium,¹ uranium,^{2,3} and plutonium.⁴ The present method has been developed for determining microgram amounts of plutonium using the Los Alamos automated spectrophotometer⁵ or a similar instrument in which all operations, from addition of reagents through the absorbance measurements, are done in a single tube. Extraction of the Pu(IV)-CPA complex into a water-immiscible alcohol improves selectivity and simultaneously increases the complex absorbance. A two-wavelength absorbance measurement improves precision by correcting for base-line shifts in a single-beam instrument, such as the Los Alamos automated spectrophotometer. The method as presented, however, provides a selective determination of plutonium suitable for a manually operated spectrophotometer.

II. EXPERIMENTAL

A. Spectra in Hydrochloric Acid and *n*-Pentanol

In dilute HCl solution, plutonium(IV) and CPA form a green complex having absorption peaks at 630 and 685 nm. The peaks of the extracted complex in *n*-pentanol are shifted to shorter wavelengths of 625 and 673 nm with a four-fold increase in absorbance (Fig. 1), similar to the effect reported by Yamamoto for the U(IV)-CPA complex extracted into 3-methyl-1-butanol.⁶

Various pairs of wavelengths, one wavelength representing the peak and the other the valley, were evaluated for absorbance measurement. The pair giving best precision was 680 and 710 nm.

B. Effect of Solvent

A reliable extraction system requires a solvent that has low viscosity, low volatility, and effective phase dis-

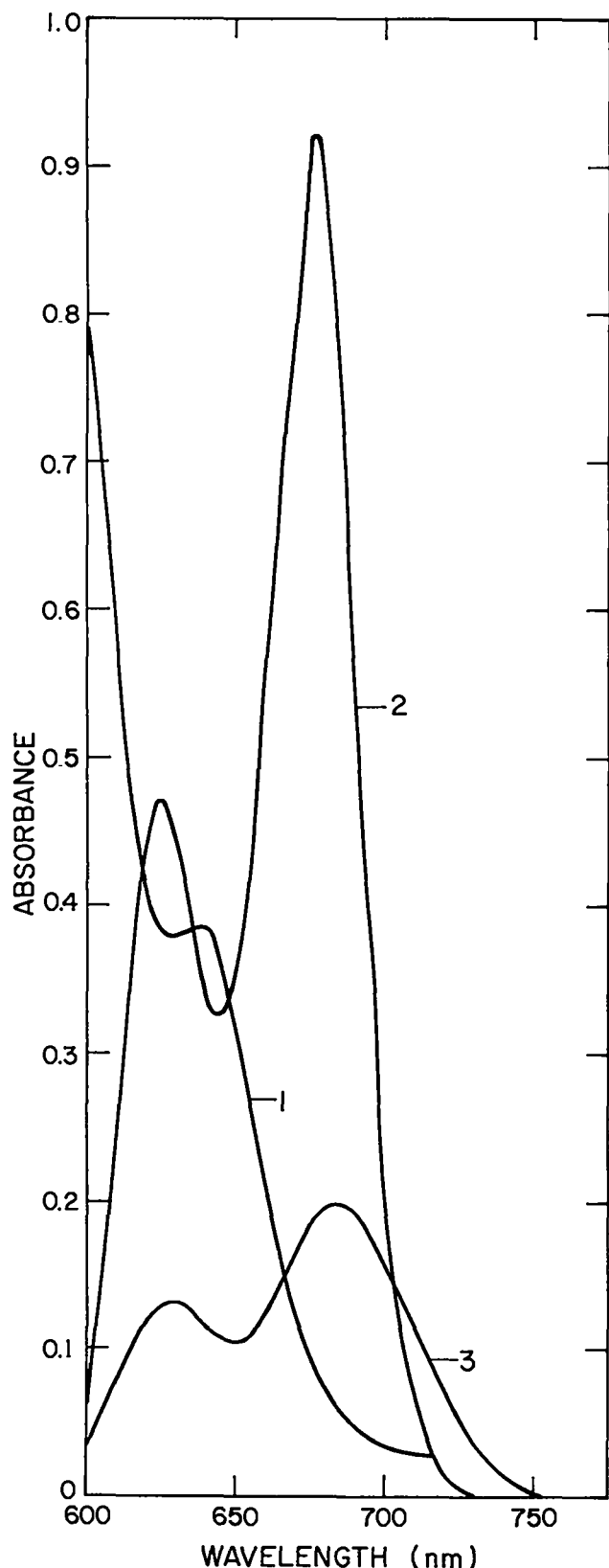


Fig. 1. Absorption Spectra. (1) 0.04% CPA in n-pentanol relative to air. (2) 15 µg Pu(IV)-CPA in n-pentanol relative to reagent blank in n-pentanol. (3) 15 µg Pu(IV)-CPA in 1.5M HCl relative to reagent blank in 1.5M HCl.

engagement. Of five water-immiscible alcohols evaluated, color intensity was higher and phase disengagement better for n-pentanol, isopentanol, and hexanol than for cyclohexanol and p-methylhexanol. N-pentanol best met the stated requirements and was selected as the extractant.

C. Effect of Hydrochloric Acid Concentration

The absorbance of the Pu(IV)-CPA complex in n-pentanol is essentially constant when extracted from 1 to 2M HCl. The solubility of hydrochloric acid in n-pentanol increases at higher acid concentrations to lower the absorbance. A concentration of 1.5M HCl was selected.

D. Effect of Plutonium Oxidation State

The extraction of plutonium, initially in the (III), (IV), and (VI) oxidation states was measured under the conditions of the recommended procedure. The Pu(III) had been prepared by drying the HI-HCl eluate from an anion exchange column, Pu(IV) by oxidizing Pu(III) with sodium nitrite, and Pu(VI) by evaporating Pu(III) to dryness with HClO₄ and a trace of HNO₃. Extraction was measured in each case by liquid scintillation counting the aqueous phase. Extractions of Pu(III), Pu(IV), and Pu(VI) were 98.0, 99.6, and 99.4%, respectively, indicating conversion to Pu(IV), which is strongly complexed by CPA. To ensure constant extraction, sodium nitrite is added in the recommended procedure to oxidize Pu(III) to Pu(IV).

E. Effect of Chlorophosphonazo III Concentration

The effect of CPA concentration on the formation of the complex was established by extracting 10-µg (0.042 µmole) portions of Pu(IV) with increasing quantities of 0.04% ($5.3 \times 10^{-4}M$) CPA solution from 1.5M HCl into n-pentanol. The absorbance, relative to a reagent blank, increased to 0.6 mℓ (0.32 µmole) of CPA and was constant thereafter. A volume of 1 mℓ of 0.04% (0.53 µmole) CPA solution was selected.

F. Effect of Time

The development and extraction of the Pu(IV)-CPA complex is not instantaneous. After 10 min of phase

mixing and a 2-min period of phase disengagement, the absorbance is 93% of maximum, increasing to 97% in 2 h and to the maximum in 12 h. This is not critical for an automated instrument operation in which the timing of all steps is highly reproducible. For manual operation, however, a 2-h color development time is recommended.

G. Measurement Reliability

Measurement precision and linearity of absorbance were determined for the range of 0 to 20 μg of plutonium. As shown in Table I, the response is linear through 17.5 μg and the relative standard deviation is better than 1.5% for the range of 2.5 to 17.5 μg .

Measurement precision of the method was obtained over a 3-month period. The relative standard deviation for 60 analyses of 10 μg of plutonium reference material was 1.7%.

H. Diverse Ion Tolerances

Diverse ion tolerances were established by adding individual ions to 10 μg of plutonium and comparing the absorbances to that of plutonium alone. The initial molar ratio tested was 1000. If a result was significantly different at the 95% confidence level relative to plutonium alone, lower amounts were tested until there was no significant difference. The results, summarized in Table II, show high tolerances for many impurities present in plutonium materials characteristic of the nuclear fuel cycle.

Additional tolerance to diverse ions is difficult to obtain with masking agents because they also form highly stable complexes with Pu(IV).⁷ Additional selectivity is gained by an ion exchange separation preceding the extraction, described in Sec. III, Recommended Procedure. The tolerances for selected ions, following the anion exchange separation, are given in Table III.

The anion exchange separation does not lower precision significantly; however, plutonium recovery decreases about 2%. Occasionally a set of samples and reference materials yield low recovery for all members of the group. This bias is caused by a residue of fine particles or degradation products of the resin that

accompanies the eluted plutonium. Wet sieving of the resin to remove fine particles usually eliminates the interference; however, it occasionally recurs. This residue resists oxidative chemical treatments including fuming with perchloric and nitric acids. It can be decomposed by ignition, a treatment not recommended for routine plutonium analyses.

III. RECOMMENDED PROCEDURE

A. Reagents

- Hydrochloric Acid (saturated, 6M, 1.5M). Distill ~6M HCl in a quartz still. Using polyethylene apparatus, pass HCl gas into the distilled acid while cooling in an ice bath. Verify that the molarity of the gas-saturated hydrochloric acid is at least 12M by titrating with a standard NaOH solution. Dilute as necessary with distilled water to obtain 6M and 1.5M HCl.
- Hydriodic Acid-Hydrochloric Acid Mixture (0.1M HI-12M HCl). Dilute 1 mℓ distilled 7.4M HI to 74 mℓ with saturated HCl. Prepare this reagent fresh before use. The HI is distilled, preferably in quartz apparatus, to remove H_3PO_2 , the holding reductant in AR grade HI. Store the distilled HI in flame-sealed ampoules or in plastic containers under inert atmosphere.
- Nitric Acid (15.6M). Use quartz-distilled HNO_3 .
- Perchloric Acid (6M). Mix 50 mℓ 12M HClO_4 with 50 mℓ distilled water.
- Anion Exchange Resin, Bio-Rad AG MP-1, 50-100 mesh, chloride form. Separate the fines by backwashing with distilled water.
- Chlorophosphonazo III (0.04% solution in 0.5M HCl). ICN Pharmaceuticals, Inc., Plainview, NY 11803. Dissolve 40 mg CPA in a 100-mℓ volumetric flask using 96 mℓ distilled water and 4 mℓ 12M HCl. This solution is stable in the absence of strong oxidizing or reducing agents.
- n-Pentanol. Kodak Laboratory Chemicals, Rochester, NY 14650. Use without further purification.
- Sodium Nitrite (1M). Dissolve 6.9 g of the reagent in 100 mℓ distilled water.

TABLE I. Linearity of Absorbance and Measurement Precision^a

Plutonium (μg)	Absorbance	Absorbance Relative to Reagent Blank per μg Plutonium	Rel. Std. Dev. (%)
0.	0.037	0.000	11.4
1.	0.091	0.054	2.8
2.5	0.170	0.053	0.7
5.	0.310	0.055	1.5
7.5	0.439	0.054	1.0
10.	0.579	0.054	0.7
12.5	0.709	0.054	1.4
15.	0.842	0.054	0.8
17.5	0.962	0.053	0.5
20.	1.056	0.051	2.0

^aBased on 6 determinations at each level.

TABLE II. Tolerances of Metals and Nonmetals without Anion Exchange Separation (Mole Ratio, Relative to 10 μg Pu)

	1000:1	100:1	50:1	10:1	1:1
Cd II	Pt IV	Cu II	Fe III	Sb III	Cr III
Co II	Re IV	Mn II		Mn VII	Cr VI
Hg II	Se IV	Mo VI			
Pd II		Ru III			
acetate	iodate	V V			
borate	nitrate	Zn II			
bromide	nitrite	sulfide			
bromate	perchlorate	thiocyanate			
chloride	peroxide				
citrate	sulfamate				
EDTA	sulfate				
fluoride	sulfite				
formate	tartrate				
iodide	thiosulfate				

Ce III, Ir III, Rh III, phosphate not tolerated at molar ratio of 1.

TABLE III. Tolerances of Metals with Anion Exchange Separation (Mole Ratio, Relative to 10 μg Pu)

100:1		10:1
Ce III	Fe III	Rh III
Cr III	Mn VII	
Cr VI	Sb III	
Ir III	U(VI)	

Phosphate not tolerated at molar ratio of 1.

B. Apparatus

- A double-beam recording spectrophotometer with 1-cm cells was used for absorbance measurements.
- An inversion rotator, 48 rpm, Scientific Industries, Inc., Springfield, MA 01103, was used for extracting the Pu(IV)-CPA complex into pentanol.

C. Ion Exchange Separation (optional)

1. Column Preparation. Prepare the columns from disposable 1.6 ml polyethylene droppers having 2 ml reservoirs (Cole-Parmer Instrument Co., Chicago, IL 60648, No. C-6096-00 or equivalent). Cut off the top of the dropper to form an open reservoir and insert a small wad of quartz wool in the tip. Transfer a water slurry of the ion exchange resin into the column until the settled height is just below the reservoir. Immediately before use, pass 6 ml of the specially prepared, saturated HCl through the column. Keep the column wet with saturated HCl until use.

2. Column Separation.

(a) To an aliquot containing between 1 and 17.5 μg of plutonium in a 10-ml Teflon[®] beaker, add 2 drops of 6M HClO₄ and 1 drop of concentrated HNO₃, then evaporate to dryness.

- All Pu(III) oxidizes.

(b) Dissolve and transfer the residue with a polyethylene dropper to the anion exchange column with three 0.5-ml portions of saturated HCl.

- Americium, thorium, and many interfering elements elute.

(c) Pass 8 ml saturated HCl through the column in

2-ml increments.

- This completes the elution of many interfering elements.

(d) Place a 15-ml glass test tube under the column and elute the plutonium fraction with two 2-ml portions of 0.1M HCl mixture, then wait 10 min.

- Plutonium reduces to Pu(III).

(e) Pass four 2-ml portions of the 0.1M HI-12M HCl mixture through the column to complete the reduction and quantitative elution of plutonium.

(f) After allowing the saturated HCl eluate from step (e) to outgas overnight, gently evaporate the solution to dryness.

(g) Add one drop of concentrated HNO₃ to the residue and heat to remove the residual HI.

(h) Convert the residue to a chloride salt by heating to dryness with one drop of 6M HCl.

(i) Redissolve the salt in 4 ml of 1.5M HCl.

D. Extraction and Absorbance Measurement

- Using the recommended procedure and standard reference material, establish a calibration curve in the range of 0 to 17.5 μg of plutonium.

1. Deliver a 4-ml aliquot containing between 1 and 17.5 μg of plutonium in 1.5M HCl into a glass test tube.

2. Add 0.1 ml 1M NaNO₂ and warm the solution in a 60°C water bath for 10 min.

- Pu(III) oxidizes to Pu(IV).

3. After cooling, add 1 ml 0.04% CPA solution and briefly agitate the mixture.

- Pu(IV)-CPA complex forms.

4. Add 10 ml n-pentanol. Seal the test tube with a hollow polyethylene stopper using masking tape and mix the two phases for 10 min using a 48-rpm inversion rotator.

- Pu(IV)-CPA complex extracts into n-pentanol.

5. After phase disengagement, measure absorbance at 680 and 710 nm. For manual operation, a 2-h wait is recommended before absorbance measurement.

- The difference between the two absorbance measurements is proportional to the plutonium concentration.

ACKNOWLEDGMENTS

Technical guidance and editorial suggestions of J. E. Rein are appreciated.

REFERENCES

1. T. Yamamoto, "Extraction-Photometric Determination of Thorium with Chlorophosphonazo III," *Anal. Chim. Acta* **63**, 65-70 (1973).
2. V. F. Luk'yanov, A. M. Lukin, E. P. Duderova, and T. E. Barabanova, "Photometric Determination of the Uranyl Ion by Means of Phosphonazo III and Chlorophosphonazo III in Strongly Acid Solutions," *Zh. Anal. Khim.* **26**, No. 4, 772-776 (April 1971).
3. T. T. Bykhovtsova and I. A. Tserkovnitskaya, "Photometric Determination of Uranium VI by Means of Chlorophosphonazo III in Organic Solutions after the Extraction Separation of Uranium by Means of Trialkylamines," *Zh. Anal. Khim.* **32**, No. 4, 745-750 (April 1977).
4. T. Yamamoto, "Spectrophotometric Determination of Plutonium with Chlorophosphonazo III," *Mikrochim. Acta* 871-877 (1974).
5. D. D. Jackson, D. J. Hodgkins, R. M. Hollen, J. E. Rein, "Automated Spectrophotometer for Plutonium and Uranium Determination," Los Alamos Scientific Laboratory report LA-6091 (February 1976).
6. T. Yamamoto, "Extraction-Photometric Determination of Uranium IV with Chlorophosphonazo III," *Anal. Chim. Acta* **65**, 329-334 (1973).
7. F. W. E. Strelow and T. N. Van der Wart, "Comparative Study of the Masking Effect of Various Complexans in the Spectrophotometric Determination of Uranium with Arsenazo III and Chlorophosphonazo III," *Talanta* **26**, 537-542 (1979).

Printed in the United States of America
Available from
National Technical Information Service
US Department of Commerce
5285 Port Royal Road
Springfield, VA 22161
Microfiche \$3.50 (A01)

Page Range	Domestic Price	NTIS Price Code	Page Range	Domestic Price	NTIS Price Code	Page Range	Domestic Price	NTIS Price Code	Page Range	Domestic Price	NTIS Price Code
001-025	\$ 5.00	A02	151-175	\$11.00	A08	301-325	\$17.00	A14	451-475	\$23.00	A20
026-050	6.00	A03	176-200	12.00	A09	326-350	18.00	A15	476-500	24.00	A21
051-075	7.00	A04	201-225	13.00	A10	351-375	19.00	A16	501-525	25.00	A22
076-100	8.00	A05	226-250	14.00	A11	376-400	20.00	A17	526-550	26.00	A23
101-125	9.00	A06	251-275	15.00	A12	401-425	21.00	A18	551-575	27.00	A24
126-150	10.00	A07	276-300	16.00	A13	426-450	22.00	A19	576-600	28.00	A25
									601-up	†	A99

†Add \$1.00 for each additional 25-page increment or portion thereof from 601 pages up.